

Sprint GT

Model 2.593

The MEYRA **SPRINT GT** **electronic wheelchair** meets mobility and enhanced independence requirements because of the MEYRA accredited technology applied to this new model. It has the benefits of the SPRINT GT, both for the home and outside.



CONTROL BOX

The LCD-screen drive control box contains a counter of the distance covered daily, a speed indicator and a 3-digit LCD combination lock as a programmable block to starting.



OBSTACLES

The front wheels are equipped with an axle system which flexes a few degrees to cope easily with any obstacle.



SEAT UNIT

The seat unit may be folded backwards for maintenance and preservation work.



FOOTREST

Easy change from /in the wheelchair. The foot area is free and there is no danger of collision.



SIDES

The side parts generally come with cushioned armrests and may be blocked and disassembled as well as adjusted vertically.



NON-STANDARD ACCESSORIES

REMOVABLE CONTROL PANEL

Removable control box for access to work table.




SPECIAL ADJUSTABLE BACK

The back is adjustable and its angle may be varied.



Dimensions and Weight

Total length with footrest and anti-roll wheels	113 cm.	Seat height	54 cm.
Length without footrests	85 cm.	Seat depth	39-50 cm.
Total width	65-67 cm.	Seat armrest	20-25 cm.
Height without arms and without back	40-48 cm.	Maximum user weight	120 kg.
Total height	102-56 cm.	Maximum permitted weight	245 kg.
Seat backrest	45 cm.	Acid batteries	

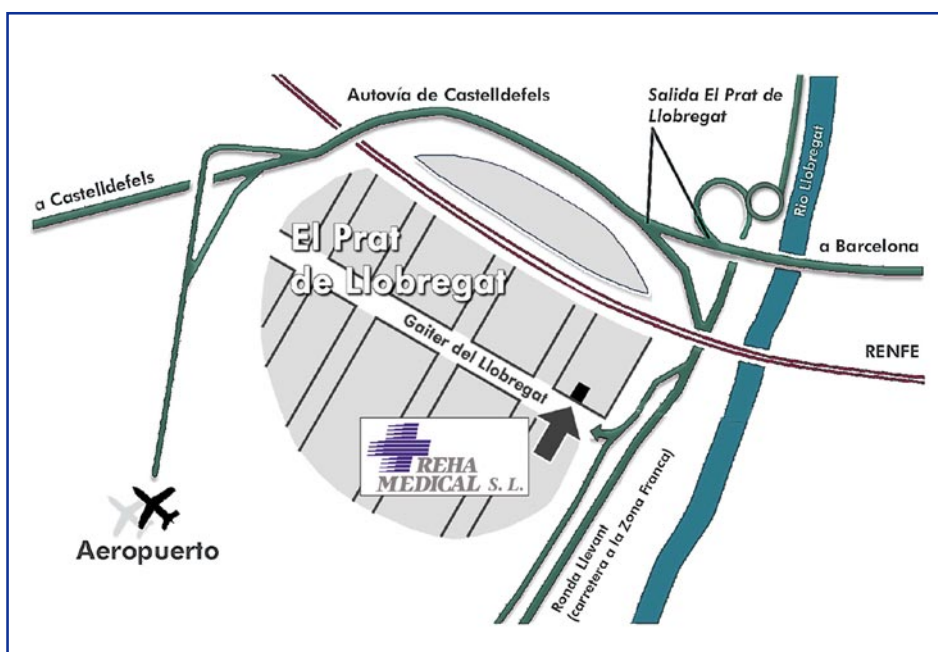

Polígono Industrial Estruc • Gaier del Llobregat, 131

Tel. / Fax 93 478 71 08

08820 El Prat del Llobregat • Barcelona (Spain)

Apdo. de correos nº 16

www.lifante.net

E-mail: reha-medical@lifante.net