

Lucca

Model 1.825

Lucca wheelchairs for children provide mobility and independence and consequent increased quality of life and enhanced safety.

Suitable for use both inside and for open air leisure activities.



FOOTRESTS

The foot area is free and there is therefore no danger of collision. It may be driven in short movements (foot-driven forward wheel chair drive).



BRAKES

According to the version, **pressure brakes** may exert direct pressure on the tyres or it may be additionally equipped with **drum brakes**. These are operated by the accompanying person and recognisable by the manual brake lever and the larger wheel hub size.



HAND GUARD / SPOKES

The hand guard / spokes protect the child's hands and prevent them from getting caught accidentally in the spokes when the chair is moving. They also prevent damage to the wheel spokes caused by impact or by small bumps.



FOLDING/UNFOLDING

The wheelchair may be folded without the use of tools in very few operations. It is easily folded and transported by car. The wheels may also be easily removed from their axle by using an automatic anchorage button.



ARMRESTS

The side pieces may be disassembled and come equipped with cushioned armrests and an elegant plastic side panel.



ANTI-ROLL WHEELS


It has an anti-roll system with two small wheels to enhance stability.



Lucca

Dimensions and Weight

Total length with footrests	77 cm.	Seat height variable	41-46 cm.
Total width	49-51-53-55-57 cm.	Seat depth	30-32-34-36-38 cm.
Seat width	28-30-32-34-36 cm.	Maximum user weight	75 kg.
Total height	70 cm.	Maximum permitted weight	88 kg.

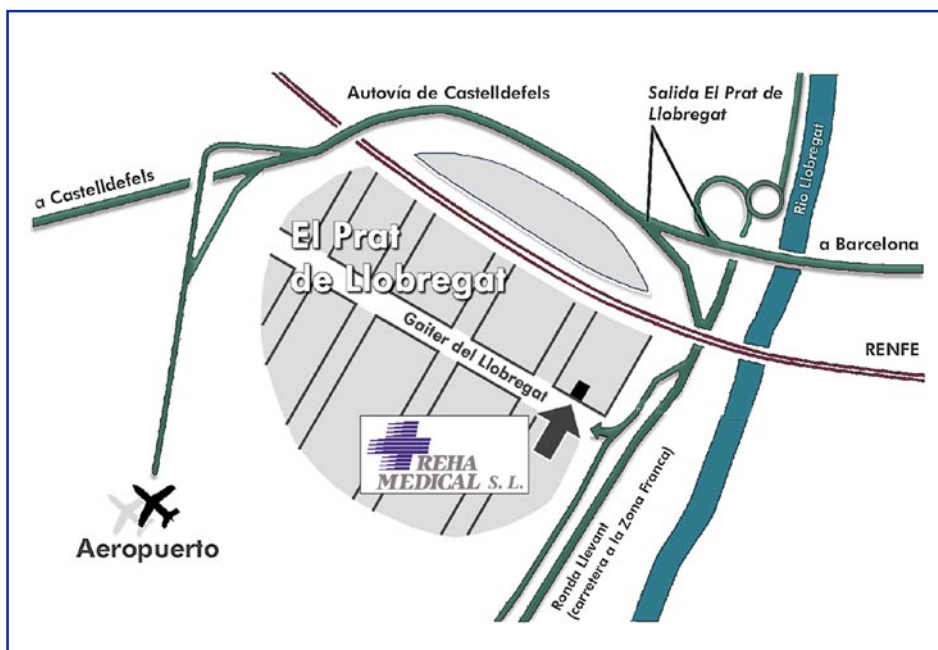

Polígono Industrial Estruc • Gaijer del Llobregat, 131

Tel. / Fax 93 478 71 08

08820 El Prat del Llobregat • Barcelona (Spain)

Apdo. de correos nº 16

www.lifante.net

E-mail: reha-medical@lifante.net